



DISINFECTION CABINETS FOR KNIVES

MZ05 - 05S - 05P - 05SP

MZ10 - 10S - 10P - 10SP

MZ20 - 20S - 20P - 20SP

TECHNICAL INSTRUCTIONS



MZ10SP



MZ20S

Thank you for buying this professional disinfection cabinet SOFINOR.

In order to obtain optimum performance from this equipment, please read the information in this brochure carefully and follow the recommendations concerning:

- Installation,
- Use,
- Maintenance.

This appliance is for professional use only and must be used by qualified personnel.

It must be set up in compliance with the regulations in force.

DECLARATION OF CONFORMITY


■ This electric appliance is made in compliance with the European Directives

- « Low Tension » **Ref. 73/23/CEE**
- « Electromagnetic Compatibility » **Ref. 89/336/CEE**

■ The applied harmonised norms are as follows :

- Security of electro-domestic appliances and analogous **EN 60 335-1**
- CEM – norm generic emission (Residential, commercial and light industry) **EN 50 081-1**
- CEM – norm generic immunity (Residential, commercial and light industry) **EN 50 081-1**

Symbol for separate collection applicable in European countries

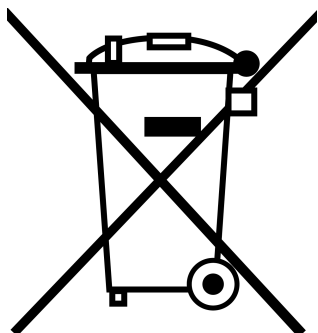

This symbol indicates that this product is to be collected separately.
The following apply only to users in European countries.

- * This product is designated for separate collection at an appropriate collection point. Do not dispose of as household waste.
- * For more information, contact the retailer or the local authorities in charge of waste management.



FICHE TECHNIQUE

MODELS :	Cabinet with magnetic rail			Cabinet with basket rail		
	5 knives 414x150x646	15 knives 580x150x646	30 knives 1060x150x646	15 knives 414x150x646	25 knives 580x150x646	50 knives 1060x150x646
Model without key	MZ05	MZ10	MZ20	MZ05P	MZ10P	MZ20P
Model with key	MZ05S	MZ10S	MZ20S	MZ05SP	MZ10SP	MZ20SP

- Made of stainless steel (AISI 304L)
- A transparent plexi-glass door with a magnetic closing or key lock
- A timer (0 – 2 hours)
- A germicide lamp 08 W (MZ05), 15 W (MZ10) or 30 W (MZ20) single phase 220-230V 50Hz (60Hz for United Kingdoms)
- Sprinkle protection: IP24
- Automatically stop if the door is open

INSTALLATION

1.1) SETTING UP

- Wall fixing with 2 screws through the metal.
- Pierce two holes in the wall, countermarked with the cabinet and place in the holes two pegs and screws without tightening.
- Hang up the cabinet on the 2 screws and make sure that it is in horizontal position before tightening the screws.


1.2) ELECTRICAL CONNECTION

The appliance must be supplied with the tension 220-230V single phase, 50Hz (60Hz for U.K)

Electrical connection which does not agree with the indicated tension “220-230V single phase, 50Hz” (60Hz for U.K) cause the cancellation of the guarantee.

- If the supply cable is damaged, it must be replaced by the special supply cable reference: 6570 MZ 3075 (6570 MZ 3100ANG for U.K) available in the factory or in the after-sales service.

USE


-  Connect the disinfection cabinet to the network 220-230V single phase according to indications. Place the knives on the magnetic bar or in the basket carrier after their cleaning and drying. **Cabinet should be used carefully. The radiation UVC could cause physiological disorder. Eyes and uncovered parts of the body should be kept away from the radiation of the germicide lamp.**
- Close the glass door(s) because it is essential to protect the skin and the eyes from the ultraviolet radiation, opened door can also disturb the disinfection of the knives.
- Set the timer on 30 min. The bacterium are eliminated by the ultraviolet radiation.
- After 30 min, the tube goes off automatically.
Return the knives on the magnet to disinfect the other side during 30 min.
It is not necessary to return the knives in case of cabinet with slides.
- Set the timer on 30 min and wait for the automatically stop.
- We recommend disinfecting the knives just before use.
- If the knives are not used immediately, put them into an alcoholic solution to keep the disinfection effect.

CLEANING

- Disconnect the cabinet before cleaning.
- The timer will be protected from the infiltration
- Do not use abrasive or corrosive products (solvents, scour powder, acid detergents)
- Do not clean with the hose.
- Use an antistatic rag or humid sponge to remove dust from the plexi-glass door. (**Proscribe all alcoholic products**)
- Use only a sponge or a soft rag with a soap water to clean all chromes or visible parts.
- Use a rag or a humid sponge to clean the non-visible parts (under cabinet).

MAINTENANCE

- Maintenance operations consist in replacing of the damaged or used parts.
- We recommend using only original SOFINOR parts that you can get from your retailer.
- The after-sale service SOFINOR remains at your disposal for all technical information concerning the replacing and delivery of parts.

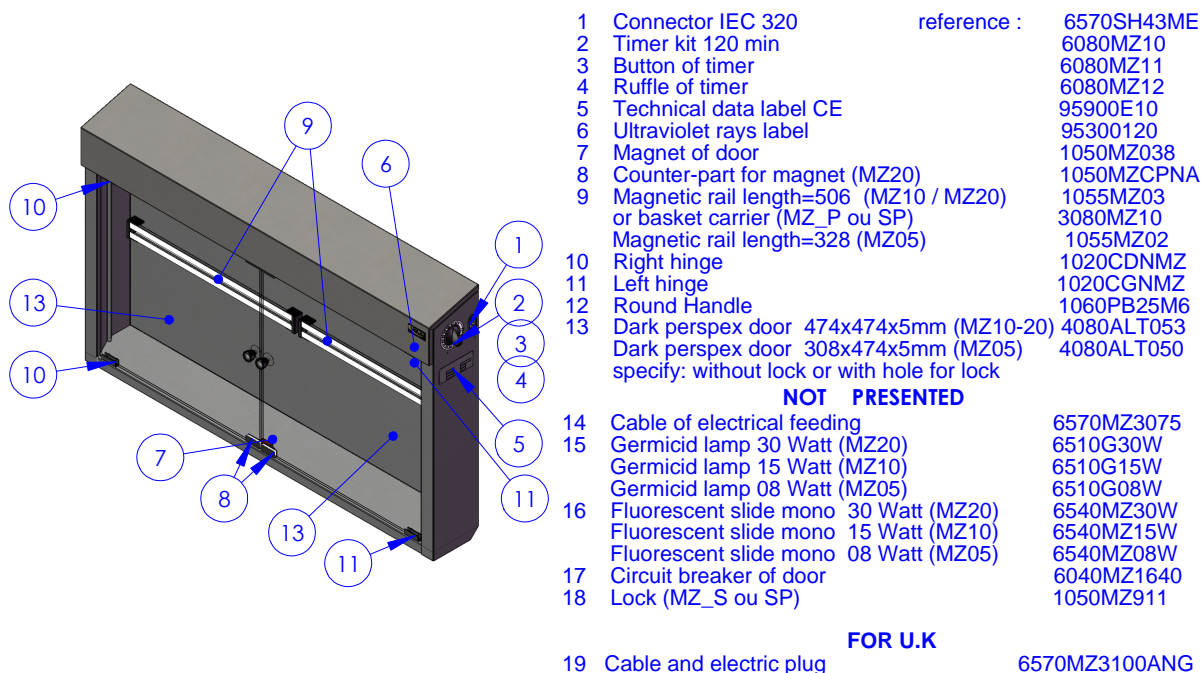

GUARANTEE LIMITS

- All parts are guaranteed 1 year.
- Our guarantee is limited to the replacing of damaged parts, in case of the normal use of the cabinet, described in this instruction.
- Guarantee is not valid for the parts, which were damaged because of:
 - an impact,
 - a connection to the wrong electric supply,
 - an use of abrasive or corrosive cleaning products,
 - a setting of non-guaranteed parts.



Z.I de la Houssoye - rue Ambroise Paré - F 59280 BOIS-GRENIER
Tel +33.3.20.10.34.10 - Fax +33.3.20.35.95.95 - www.sofinor.com
TECHNICAL INTRUCTIONS- DISINFECTION CABINETS

N 26.01.15

Page 4 / 4